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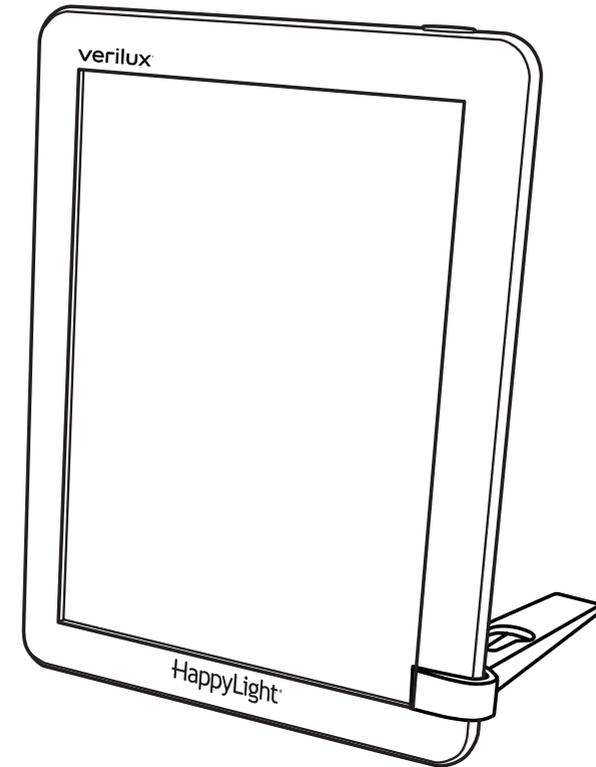
**Model: VT22**

Dear Customer,

Thank you for purchasing the Verilux® HappyLight® Lucent™ Series Energy Lamp. You now own an innovative product, manufactured to the highest standards and backed by a one-year limited warranty.

This product brings Natural Spectrum Daylight into your home or office. Many other healthy-lighting products are available online. Visit us on the web at [www.verilux.com](http://www.verilux.com) to learn more about all of our quality Verilux products, or call our toll-free number: 1-800-786-6850.

As a Verilux customer, your satisfaction means everything to us. We look forward to serving you now and in the future.

Have a bright day!

The Verilux Customer Service Team

## Contents

FCC Notice . . . . .	3
Important Safeguards . . . . .	4-5
Setup . . . . .	5-6
How to Use . . . . .	7
Positioning . . . . .	7
Device Pretest . . . . .	7
Establishing a Daily Routine . . . . .	8
Features and Benefits . . . . .	8
Care and Cleaning . . . . .	9
Troubleshooting . . . . .	9
Technical Specifications . . . . .	9
Warranty . . . . .	10

## FCC Notice

Model: HappyLight® Lucent LED Light Therapy Lamp (VT22)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Verilux  
2803 S. Yale St.  
Santa Ana, CA 92704  
[verilux.com](http://verilux.com)

The product is protected against unauthorized use. DO NOT modify it without the authorization of the manufacturer, or the user's authority to operate could be voided.

**Note:** This equipment has been tested and found to comply with the limits for a Class B<sup>†</sup> digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

<sup>†</sup> This Class B digital apparatus complies with Canadian ICES-005.

## READ ALL INSTRUCTIONS BEFORE USING

### Important Safeguards

#### ⚠ DANGER:

- To avoid electrocution, do not operate this device near water.

#### ⚠ WARNING:

- Close supervision is recommended when these lamps are used by or near children, or others who may not understand the need for these precautions.
- To prevent risk of shock or personal injury when cleaning lamp, make sure you have turned off and unplugged the lamp and allowed sufficient time for the lamp to cool.
- Do not shorten or cut the power cord.
- Certain medications (e.g., anti-depressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity or hypomania. Consult with your healthcare provider before using this product.
- If you have a pre-existing eye condition or are at risk for developing Age Related Macular Degeneration or other eye conditions, consult with your healthcare provider before using this product.
- Some people experience light sensitivity to bright light or upon over-exposure to bright light. Discontinue use if you experience any discomfort.

#### ⚠ CAUTION:

- This lamp is for indoor use on a hard surface (e.g., a table or desk) only.
- Do not operate this lamp if it has been damaged in any way.
- Do not disassemble. There are no user-serviceable parts in this lamp.
- Use only the power adapter supplied with your lamp. If other power plugs, cords or connectors are used, damage to your lamp or power adapter may occur.
- This lamp may cause interference with radios, cordless telephones or devices that use a wireless remote control, such as televisions. If interference occurs, move the lamp away from the device, plug the lamp or the device into a different outlet or move the lamp out of the line of sight of the remote control receiver.†
- This lamp is not for use with light dimmers, timers, motion detectors or voltage transformers.

† This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.”

This Class B digital apparatus complies with Canadian ICES-005.

### Important Safeguards — continued


#### ⚠ CAUTION: — continued

- Use a quality surge suppressor on all electronic equipment. Voltage variations and spikes can damage electronic components in any system. A quality suppressor can eliminate the vast majority of failures attributed to surges. Suppressors may be purchased at electronics stores.
- Never cover the lamp or place anything on top of it when it is in operation.

### SAVE THESE INSTRUCTIONS

### Setup and Take Device Pretest


#### 1 Setup – Remove Stand

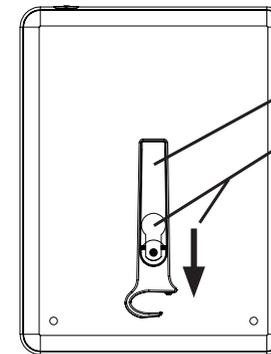- Stand affixed to back of lamp
- To remove stand – slide stand downward to remove from keyhole

#### 2 Setup – Affix Stand - vertical


Affix stand by clipping to bottom left or bottom right — as shown.

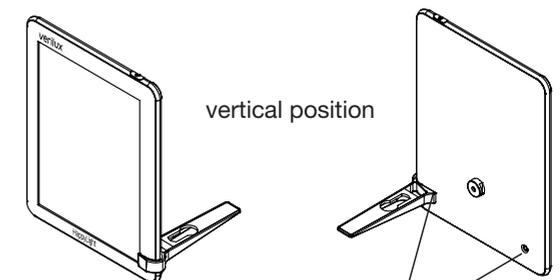Affix to indentations for secure position


### 3 Setup – Affix Stand - horizontal

Stand can be attached so lamp orientation is horizontal.

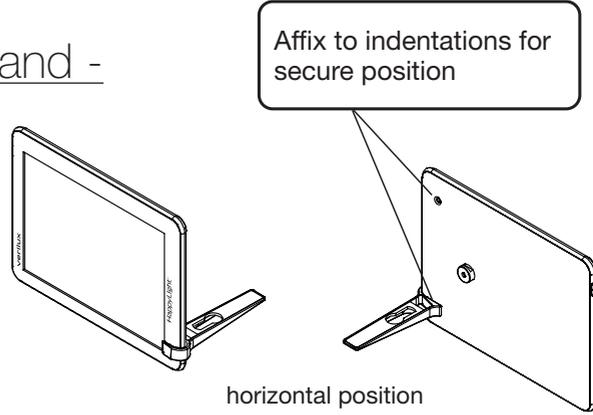horizontal position

### 4 Setup – Plug In

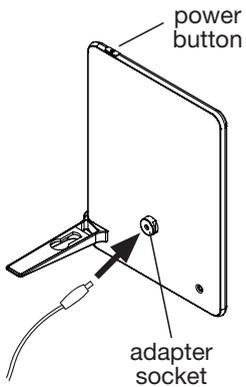- Place lamp on a flat level surface
- Insert power adapter plug into the adapter socket — as shown on the rear of the HappyLight® Therapy Lamp
- Plug adapter in to a standard 110 AC electrical outlet
- Turn on the lamp

### 5 Take Device Pretest


Before you begin using your HappyLight®, please take the Device Pretest on page 7 and read through the manual for helpful information and safety precautions.

### 6 How To Use – Positioning

Sit close to your HappyLight® Therapy Lamp with your eyes open. While staying oriented toward the light, feel free to engage in other activities such as reading, writing, hobbies and crafts, sewing, needlepoint, applying make-up or working at your computer or desk.

The lamp could be positioned as close as 6 inches but within 24 inches of your face. It is most effective when it is positioned slightly off-center. (See diagram.) *Note: Do not look directly into the lamp.*

The light from the lamp must be directed at your eyes; therefore, your eyes must be open to achieve the full benefit. *Note: Tinted glasses reduce the amount of light reaching your eyes.*
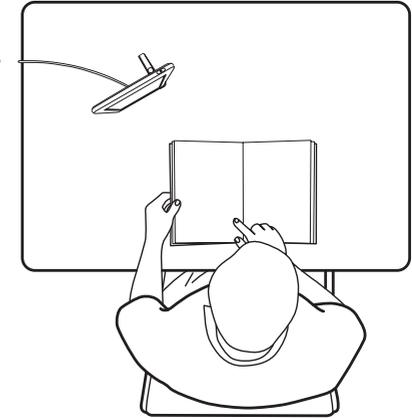
### Device Pretest

To help determine your degree of sensitivity to light, we recommend that you perform this pretest upon the first use of your HappyLight® Therapy Lamp.

Step 1 - Turn on the lamp and position it about 45cm/18 inches from your face. (Do not look directly into the lamp.)

Step 2 - Operate the lamp for approximately 10 minutes and turn it off. If the light has caused any discomfort, stop using it and consult a doctor. If you have no discomfort, you may consider gradually decreasing the distance from the lamp and increasing the use time, which results in more lamp intensity and greater effect. For your initial sessions, you should monitor the degree of any discomfort and accordingly adjust the time used and lamp location.

**⚠ WARNING:** Certain medications (e.g., anti-depressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity or hypomania. Consult with your healthcare provider before using this product.

**⚠ WARNING:** If you have a pre-existing eye condition or are at risk for developing Age Related Macular Degeneration or other eye conditions, consult with your healthcare provider before using this product.

## Establishing a Daily Routine

20 minutes  
to  
1 hour/day\*

Researchers at the National Institute of Health demonstrated that white light, in a range of 2,500 to 10,000 lux (a measurement of light intensity), helped to reverse the symptoms of “Winter Blues.”

We all have a different sensitivity to light, so it is important to find the settings, distance and time

that are comfortable for you. We recommend that you use your HappyLight® Therapy Lamp regularly and at about the same time each day. It is not necessary to sit in front of the lamp continuously. You can divide the time into several mini-sessions. Continue to use the lamp daily to feel an improvement in your mood and activity level.

When you are feeling sluggish and lethargic, sit closer to the lamp and/or spend more time with the lamp. If you feel a sort of edginess, as though you have consumed too much coffee, it may be time to turn off the lamp.

\* Recommended time is a function of brightness and distance. Individual results may vary. Also, refer to the warnings at the bottom of the Device Pretest section on the previous page.

## Features and Benefits

Your new HappyLight® Therapy Lamp is a remarkable product designed to offer unparalleled customization for maximum results and comfort. It:

- Delivers full spectrum light to your home, office, dorm room or apartment.
- Cues your body’s own energy enhancers to make you more alert and productive.
- Has up to 10,000 lux† to lift spirits, sharpen concentration and improve energy levels.
- Is designed with a compact shape and small footprint to make it easy to place in a convenient location.
- Uses Flicker Elimination Technology® and Optix® Glare Control to ensure the light is comfortable for your eyes.

†10,000 Lux (photopic measurement) at 6".

## Care and Cleaning

Your HappyLight® Therapy Lamp is made from high-quality materials that will last for many years with minimum care. You may want to periodically clean the lamp using a mild non-abrasive cleaner and soft cloth. When cleaning, make sure you have turned off and unplugged the lamp and allowed sufficient time for it to cool.

## Troubleshooting

If you are experiencing an issue with your HappyLight® Therapy Lamp, please:

- Make sure the power adapter and plug are inserted fully and securely.
- Make sure there is power to the wall outlet or try another outlet – one that is not connected to a dimmer.

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## Technical Specifications

### HappyLight Lucent LED Light Therapy Lamp

**Environment:** Operating Temperature: 10 to 35°C (50 to 95°F) Storage Temperature: -20 to 50°C (-4 to 122°F) Relative Humidity: 10 to 90% non-condensing

**Power:** Input Power to the Adapter: 100-240 VAC, 0.4A, 50/60Hz  
Input Power to the Unit: DC19.2V, 0.65A

**ETL Listed in Canada and the USA**

## Three-Year Limited Warranty

As the leading innovator of light therapy products for over two decades, Verilux, founded in 1956, backs this HappyLight® with a three-year limited warranty.

**Attention: Once opened, please do not return this product to the store where it was purchased for repair or replacement.**

Many questions can be answered by visiting [verilux.com](http://verilux.com), or you may call our Customer Happiness Team at 1-800-786-6850 during normal business hours.

This limited warranty is provided by Bear Down Brands LLC, 2803 S Yale St., Santa Ana, CA 92704.

Verilux warrants this product to be free from defects in material and workmanship for a period of three years from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux will replace defective parts of this product when available at no charge to the customer, subject to these limitations: This limited warranty does not include any postage, freight, handling, insurance or delivery fees. This warranty does not cover damage, defect or failure caused by or resulting from accident, external destruction, alteration, modification, abuse, misuse or misapplication of this product.

This warranty does not cover damage to the product resulting from return shipping or handling. Shipping insurance is recommended to help protect your product.

**Return Authorization is required for all returns.** To obtain a Return Authorization, please contact the Verilux Customer Happiness Team.

If, during the first three years of ownership, this product fails to operate properly, it should be returned as specified at [verilux.com/warranty](http://verilux.com/warranty) or as instructed by a Verilux Customer Happiness Specialist.

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

Please visit our website at: [verilux.com](http://verilux.com) or contact our Customer Happiness Team at


**1-800-786-6850**    [hello@verilux.com](mailto:hello@verilux.com)

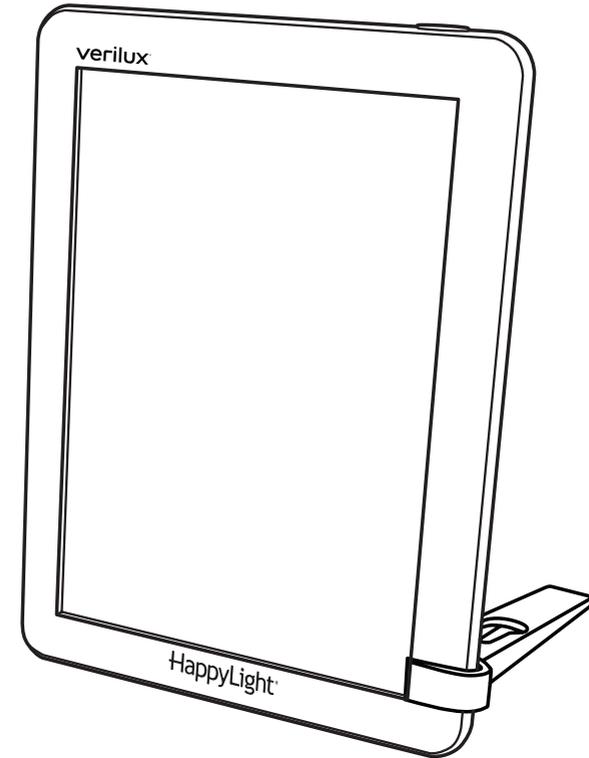Specialists are available Monday — Friday, 8am — 5pm PST

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